




NON ABSORBABLE SUTURE

SUTURES

ARTHROV[®]LINK

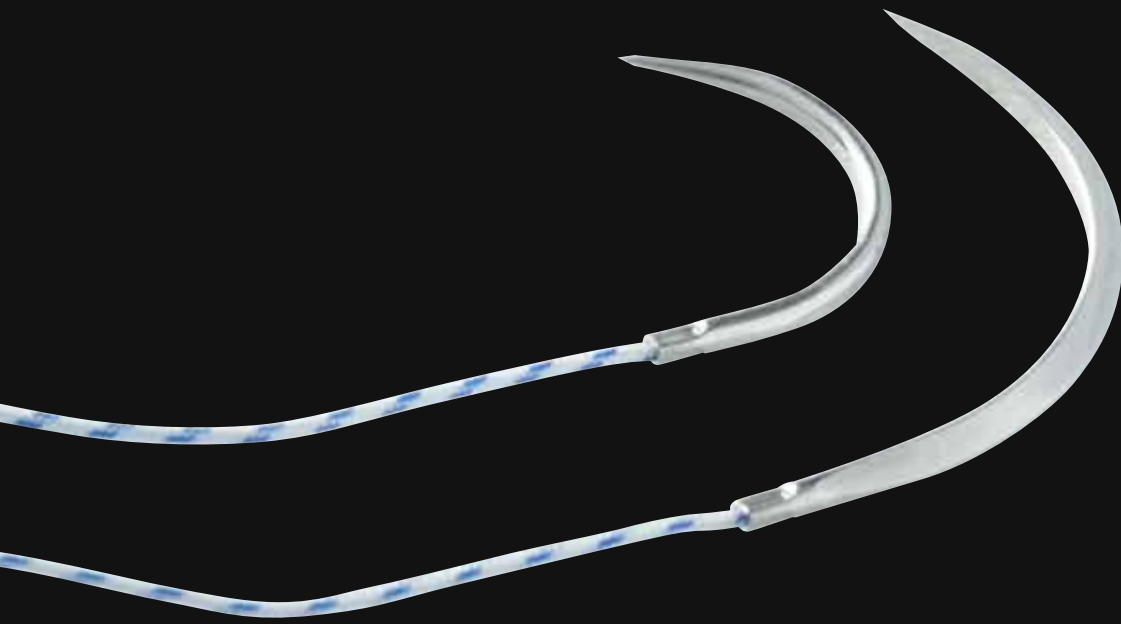

 **SPORT
MEDICINE**
ADVANCED TECHNOLOGY

NON ABSORBABLE SUTURE

For tendon and ligament reinforcement and reinsertion.

ARTHROLINK is a non-absorbable braided sterile suture prepared from threads of ultra-high molecular weight, long-chain linear polyethylene (UHMWPE). It is indicated for repair or reinforcement of ligament, tendon and ligament reinsertion, tuberosities reinsertion, approximation and/or ligation of soft tissue.



SUTURE NEEDLE Suture with double needle

Type :
Suture USP2

Size :
Length 90 cm

Colors :
white/blue – white/black – white/green – blue.

Needles :
Taper point 1/2 circle – length: 26 mm (transtendinous)
Reverse cutting 3/8 circle – length: 40 mm (transosseous)



White / Blue



White / Black



White / Green



Blue

FIBRE STRUCTURE*

Inter-locking core technology provides all fiber configurations with a sturdy core within center of the suture to act as a backbone for better knot tying and handling characteristics. This unique braid construction contributes to its superior strength. Its design is supple, soft for a better handling.





SUTURE TAPE

Tape

Type :
Suture USP2

Type :
Suture USP2

Size :
Tape length: 19 cm
Tape width: 2 mm
Suture Overall length: 100 cm

Size :
Tape length: 60 cm
Tape width: 1,6 mm
Suture Overall length: 100 cm

Colors :
blue/blue – white/black

Colors :
white/green – white/blue



SUTURE LOOP

Suture loop with needle

Type :
Suture USP2


Size :
Length 50 cm

Needles & colors:
Taper point 1/2 circle – 26 mm- white/green
Taper point 1/2 circle – 40 mm- white/blue


TENSILE STRENGTH

ARTHROLINK offers excellent mechanical properties much stronger than polyester material and market leader. The braiding allows a better handling together with a high knot strength. Sutures remain stable after in-vivo implantation.

TENSILE STRENGTH*

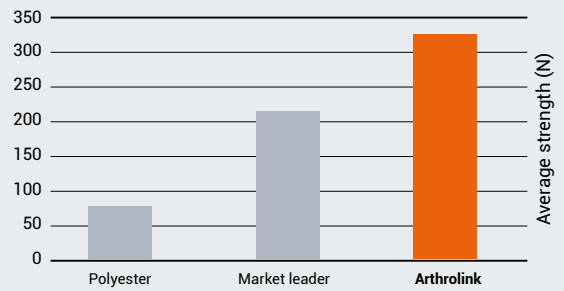

KNOT TENSILE*

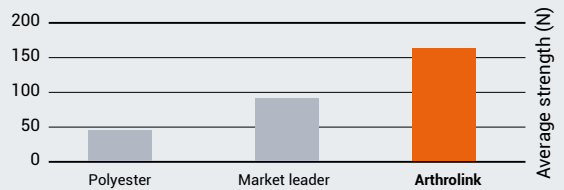

SUTURES ADVANTAGES

- HIGH TENSILE STRENGTH
- KNOT SLIDING
- SOLIDITY AND SECURITY OF THE SUTURES
- BIOCOMPATIBLE
- 4 COLORS FOR A BETTER VISIBILITY



STORAGE

ARTHROLINK must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

PRODUCTS

ARTHROLINK is packaged ready to use in a sealed pouch and sterilized by ethylene oxide. Single use. Do not re-sterilize.

For any further information, please refer to the IFU.

Designation		Colors	Reference
ARTHROLINK Suture Needle		white / blue	HSe-2C-VIMS
		white / black	HSe-2K-VIMS
		white / green	HSe-2G-VIMS
		blue	HS1008
ARTHROLINK Suture Loop		white / blue	HS955
		white / green	HS956
ARTHROLINK Suture Tape	1,6mm	white / blue	HS1007
	1,6mm	white / green	HS1006
	2,0mm	blue / blue	HS957
	2,0mm	white / black	HS958



*Data provided by the manufacturer

Document design date : Mai 2019 - Range : Non absorbable suture -
 Manufacturer : RIVERPOINT - Brand Name : ARTHROLINK - Medical Device
 Class : III - CE Mark N° : 0459 User : Health professional - Indications :
 Tendon and ligament reinforcement and reinsertion - Recommendations for
 use: it is highly recommended to read the labels and the IFU.



VIMS SA
 10 AVENUE DE FONTRÉAL - EUROCENTRE
 31620 VILLENEUVE LES BOULOC - FRANCE
 TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM